


One Machine Connection

CT Connection

MiXLink workstation is connected to the ACR/NEMA CT machine using standard network cable. The workstation can support any type of network connections, the workstation is provided with a FTP Server to accept ACR/NEMA CT images from the CT machine and send the images to dicom server using DICOM 3.0 protocol.

MRI Connection

MiXLink workstation is connected to the ACR/NEMA MRI machine using standard network cable. The workstation can support any type of network connections, the workstation is provided with a FTP Server to accept ACR/NEMA MRI images from the MRI machine and send the images to dicom server using DICOM 3.0 protocol.

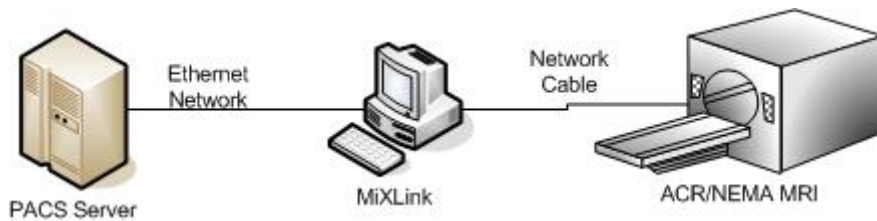

Network Connection

MiXLink workstation is connected to the ACR/NEMA MRI and CT machine at the same time using standard network cable. The workstation is provided with a FTP Server to accept both ACR/NEMA MRI and CT images from the MRI and CT machine and send the images to dicom server using DICOM 3.0 protocol.

